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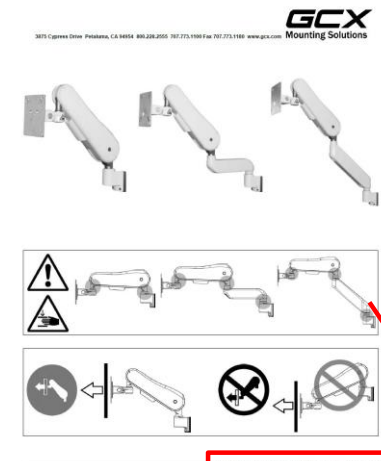
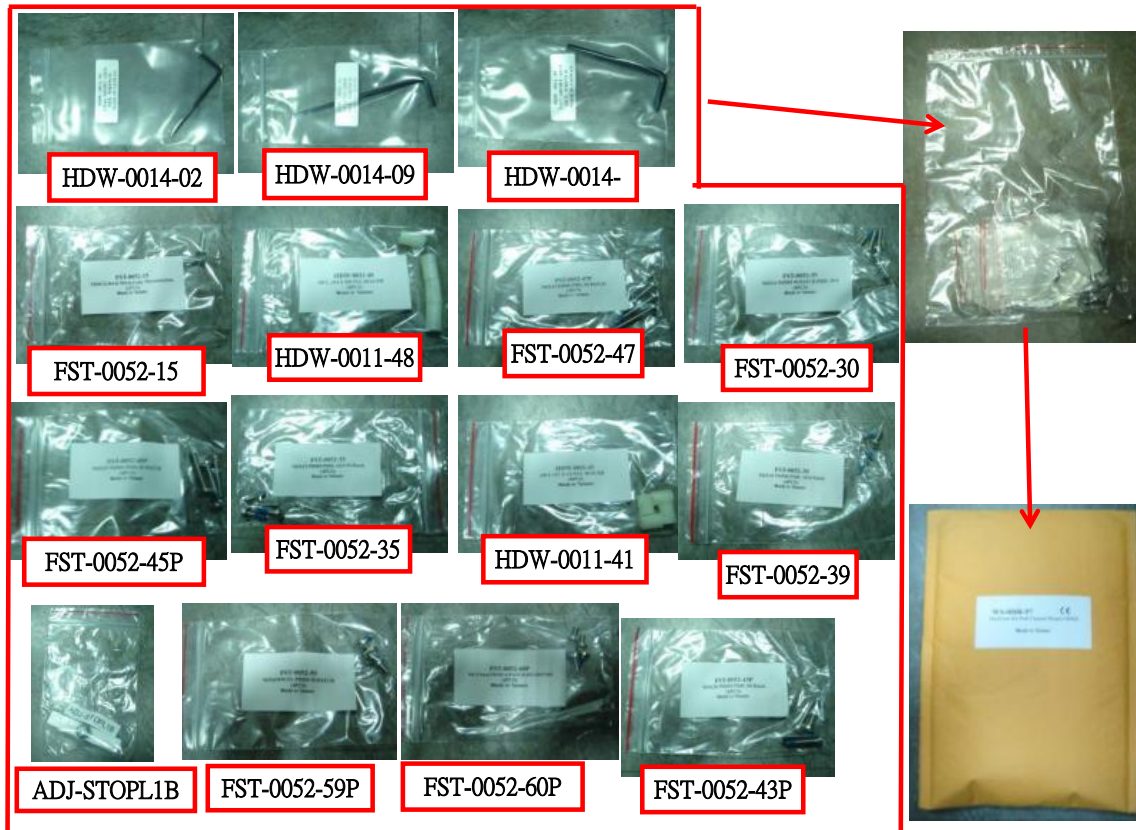
Pyramids QC Final Inspection Standard

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DU-WS-0008-101 Rev A 2/23/12 GCX Corp. Page 1 of 2

DU-WS-0008-101 REV A

Install Time: 10-15 minutes

The purpose of this manual is to describe general installation and adjustment procedures for VHM-25™ Series Arms. This manual should be used in conjunction with any instrument-specific installation guides used in mounting the monitor or display. Please read this manual and all instrument-specific installation material before installing or using this product.

### Table of Contents

1.0 Parts Reference.....	1
2.0 Tools Required.....	2
3.0 Installation and Maintenance Warnings.....	2
4.0 Installing the VHM-25 Arm in a GCX Channel.....	3
5.0 Mounting a Display on the VHM-25 Arm.....	4
6.0 Adjusting the VHM-25 Arm.....	5
7.0 Cable Management.....	7
8.0 Installing the Channel Cover.....	8
9.0 Routine Maintenance.....	9
10.0 Troubleshooting the VHM-25 Arm.....	10

### 1.0 Parts Reference

The following parts and hardware are included with this installation kit and labeled accordingly:

Item #	Description	Qty
1	VHM-25 Channel Mount Arm	1
2	1/8", 20/27, and 2.5mm Hex Wrench	1 Ea.
3	M4 x 6, 8, 10, 12, 16, 20, 25, 30mm Pan Head Machine Screws (PHMS)	4 Ea.
4	Nylon Spacers (3/8" x 1/2") and (3/8" x 5/8")	4 Ea.
5	Adjustable Stop	1
6	16 in. (40.6 cm) Wall Channel Cover (includes instructions)	1
7	M4 x 16mm Flat Head Socket Cap Screw (PHSCS)	2

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Pyramids QC Final Inspection Standard

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**GCX SPECIFICATIONS**

VHM25 IT TEST – see VHM-IT test for setting up the channel fixture.  
ATTN: The tolerances are based on the direction of the monitor as shown. LEFT or RIGHT must always tilt towards the center as where the arrows are pointed in the illustration.

LEFT	CENTER (either direction)	RIGHT
MUST ONLY READ BETWEEN $0^{\circ}$ - $2^{\circ}$	MUST BE BETWEEN $\pm 2^{\circ}$	MUST ONLY READ BETWEEN $0^{\circ}$ - $2^{\circ}$
MUST ONLY READ BETWEEN $0^{\circ}$ - $1^{\circ}$	MUST BE BETWEEN $\pm 1^{\circ}$	MUST ONLY READ BETWEEN $0^{\circ}$ - $1^{\circ}$
MUST ONLY READ BETWEEN $0^{\circ}$ - $2^{\circ}$	MUST BE BETWEEN $\pm 2^{\circ}$	MUST ONLY READ BETWEEN $0^{\circ}$ - $2^{\circ}$

PLACE DIGITAL LEVEL ON THE FRONT END YOKE AS SHOWN.

ATTN: VHM25 WITH EXTENSION – extension must be extended straight as shown.

VHM IT TEST IS THE SAME. PLEASE FOLLOW THE CHART

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DU-WS-0008-01 Rev B 9/23/13 GCX Corp. Page 1 of 10